



MANUAL

CUSTOMIZE
ACTIVITIES



Activity types

Games



Generators




Worksheets



In NeuronUP you **can customize games and generators**. Worksheets are not customizable.


1. Choose and open an activity (game or generator)

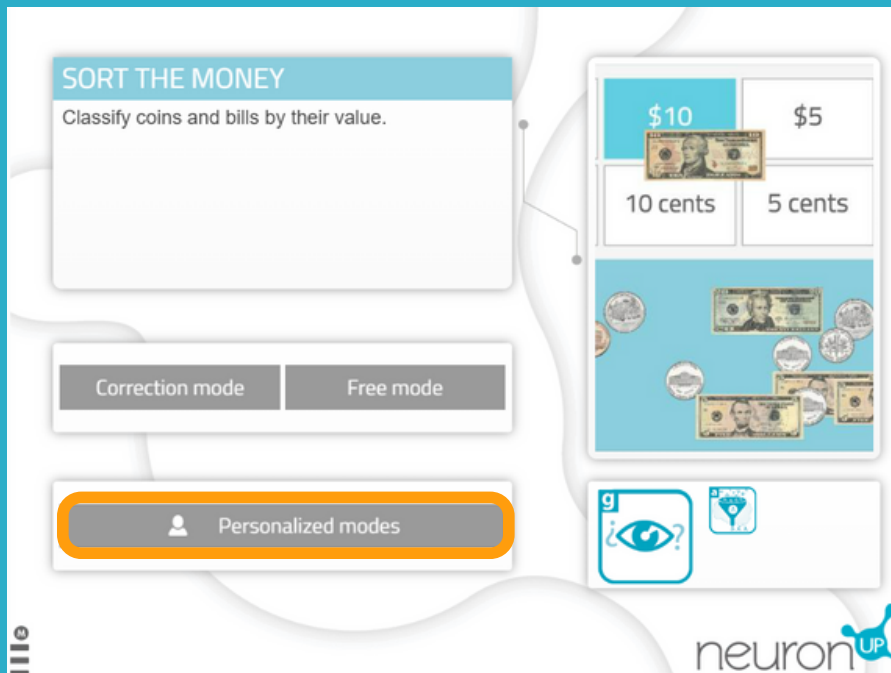

Choose a game or a generator from the NeuronUP Working Area.

Once the activity is open, select **"Personalized Modes"** for generators or **"Personalized Stages"** for games.

2. Create or edit a customization

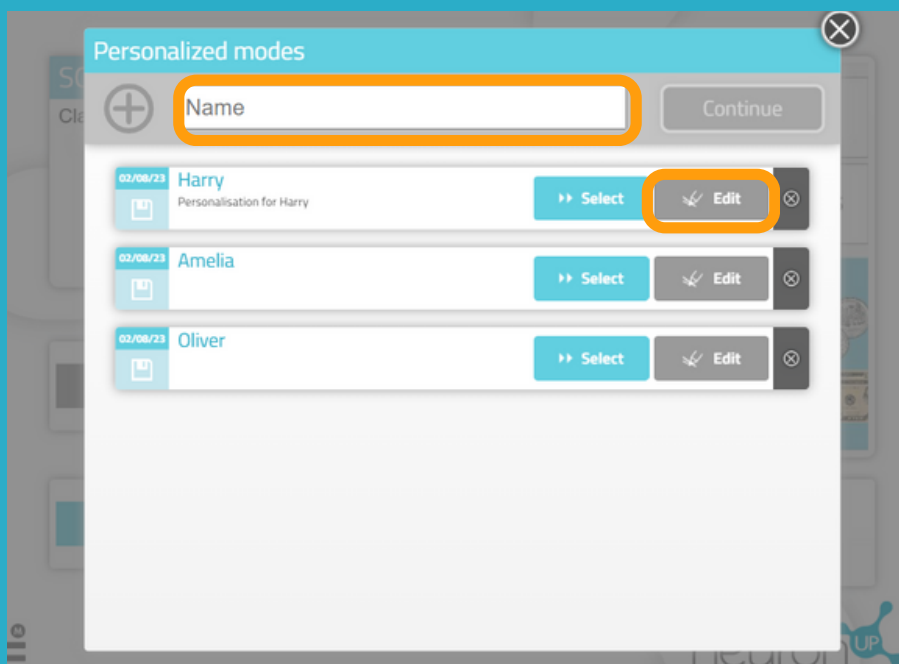

To create a new customization, **enter a name for the customization**. For example, the name of the patient it will be used for.

Then click on **"Continue"**.

To modify a customization, click on **"Edit"**.

3. Start customizing

General features

The 'General features' tab contains the following settings:

- Operation:** Radio buttons for 'Correction' (selected) and 'Free'. A text input for 'Max. errors' is set to 3.
- Maximum time:** Radio buttons for 'Yes' (selected) and 'No'. Time inputs are set to 10 Minutes and 0 Seconds.
- Inactivity warning:** Radio buttons for 'Yes' (selected) and 'No'. Time inputs are set to 2 Minutes and 0 Seconds.
- Instructions:** A text area containing the instruction: 'Put the coins and bills in their places.'
- Show timer:** Radio buttons for 'Yes, count down', 'Yes, count up', and 'No' (selected).

A 'Finish' button is located at the bottom right.

The "general features" are the same for all games and generators.

You can choose the work mode, adapt the time, warn the patient when he/she is inactive, modify the activity instructions and add a timer.

Parameters

The 'Parameters' tab contains the following settings:

- Type of currency:** A dropdown menu set to 'United States dollar'.
- No. of elements:** A text input set to 5.
- Bills/coins:** A grid of checkboxes for selecting bill and coin denominations.

Bills		
<input type="checkbox"/> 100	<input type="checkbox"/> 50	<input checked="" type="checkbox"/> 20
<input checked="" type="checkbox"/> 10	<input checked="" type="checkbox"/> 5	<input type="checkbox"/> 2
<input checked="" type="checkbox"/> 1		
Coins		
<input type="checkbox"/> 1	<input type="checkbox"/> 50 cents	<input checked="" type="checkbox"/> 25 cents
<input checked="" type="checkbox"/> 10 cents	<input checked="" type="checkbox"/> 5 cents	<input checked="" type="checkbox"/> 1 cents

A 'Finish' button is located at the bottom right.

The "parameters" are specific to each activity.

In this case, we are in a money handling activity, so you can choose the type of currency, the amount of bills/coins and their value.

Adaptability

The 'Adaptability' tab contains the following settings:

- Movement of elements:** Radio buttons for 'Drag' (selected) and 'Click'. Each option has a corresponding visual icon showing how elements are moved.

A 'Finish' button is located at the bottom right.

The "adaptability" tab allows you to choose between dragging or clicking on elements to make it easier for patients with motor impairments or those who are not good with the mouse or touch screens.

It is not available for all activities.

4. Save the customization

The screenshot shows a dialog box with the following fields and options:

- Name:** Harry
- Description:** Personalisation for Harry
- Add button:** A button with a plus icon and the text "Add".


The background interface shows tabs for "General features", "Parameters", and "Adaptability". The "Parameters" tab is active, showing options for "Operation" (Correction, Freestyle), "Maximum" (Yes, No), "Inactivity" (Yes, No), and "Seconds" (No).

You can add a **description** to the customization.

This will make it easier for you to remember its characteristics.



To save the customization, click on "Add".

The customized activity is ready and saved

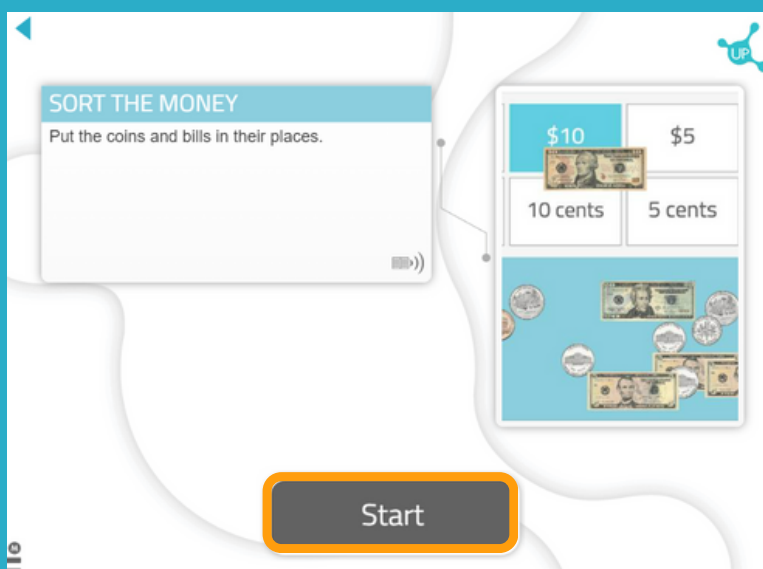


You can **work directly** with the customized activity.

When you want to work with it later, you will find it in the **list of activity customizations**.

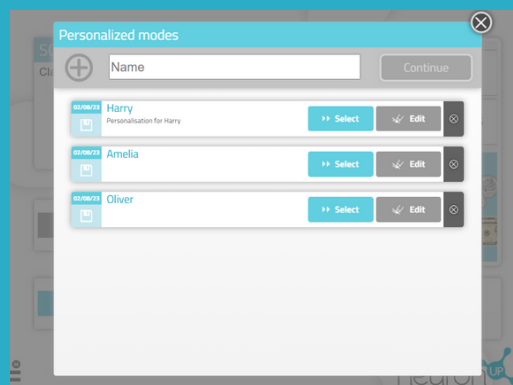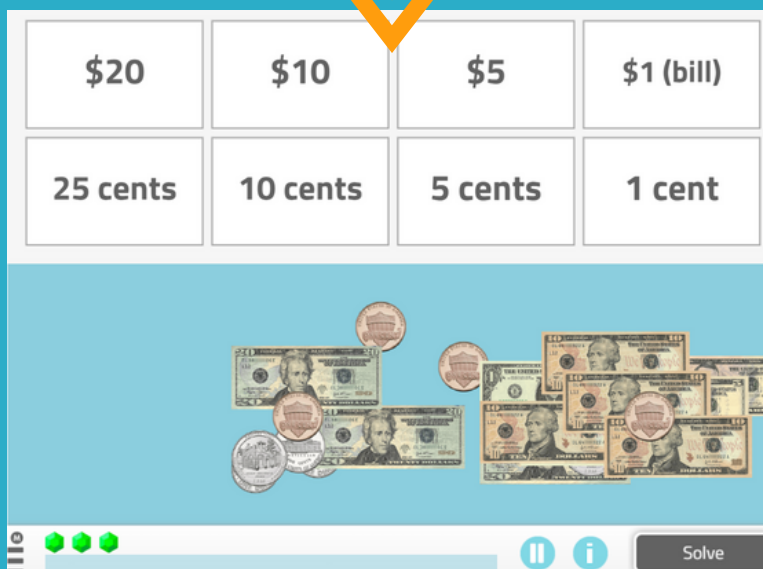

Find your customization when you create a session

1. Go to "Sessions" > "Digital Sessions" > "New Digital Session".
2. Set up your session and add the activities by clicking on them.
3. Choose the activity for which you have created a customization. In this case *Sort the Money*.

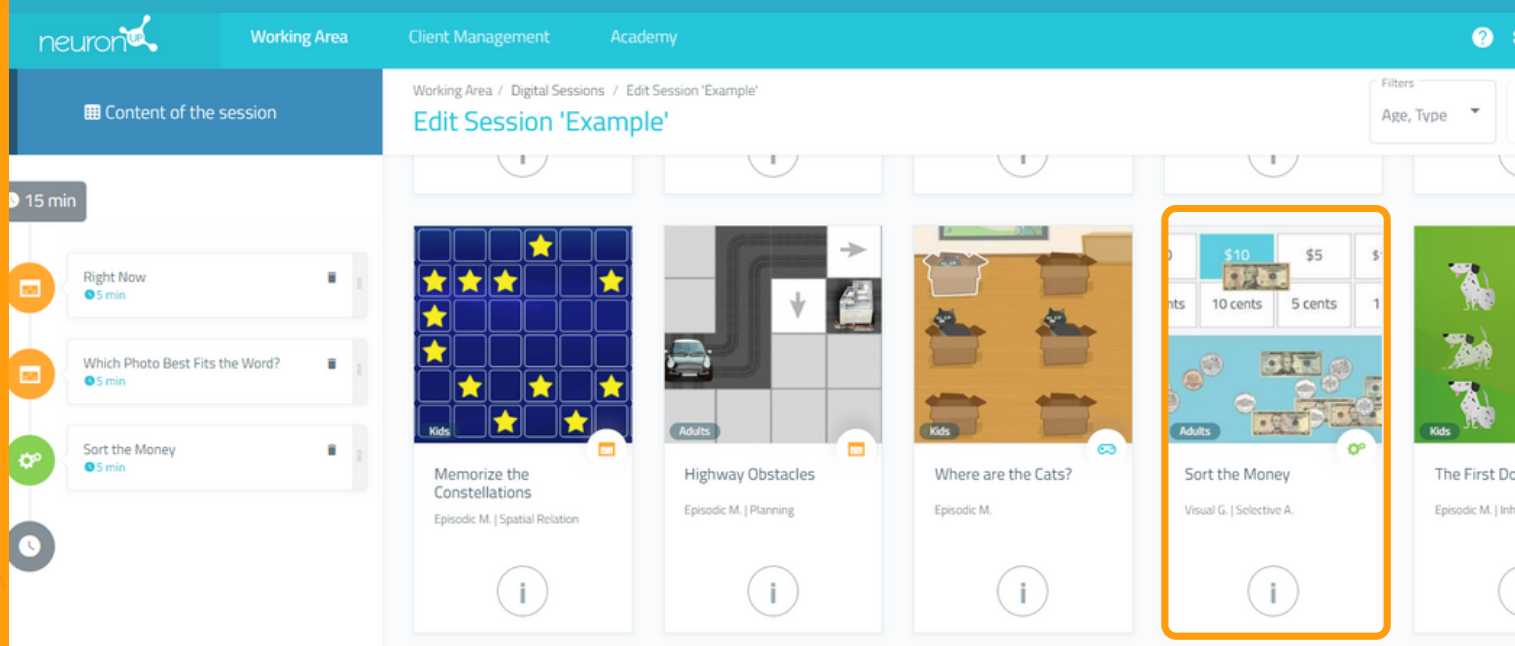

4. Once the activity has been added, click on it on the left side of the screen, as in the following image.

5. Next, activate the "Choose where to start" filter and you will see the different customizations you have created for this activity.

6. Choose the customization you want. In this example we chose "Amelia". This way, your patient will carry out the session with the customized activity.

